Light quantum

Newface
-NewFace NewLife-

NV-LD100
MANUAL BOOK

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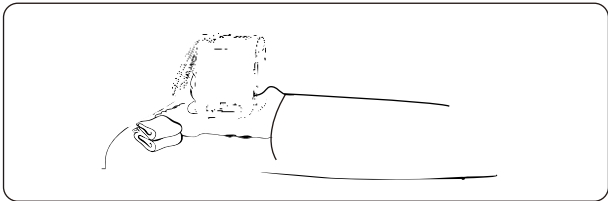
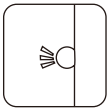
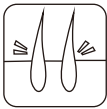
Using notice


Notice when use this machine

- Please read the product specification carefully before you use the product
- Please use the standard charging fittings
- If there is some abnormality in your skin when you use this product, please stop using and contact the doctor in time
- Please don't use this product in the bathroom and other places of high humidity
- Please don't get close to the fire
- Please don't put the product in the hot and wet environment
- Don't keep it in the confined space for a long time when it is energized
- Please don't take it apart randomly, press it or transform it
- Please don't expose it to the hot environment of above 60 degrees
- Please don't drop this product from the high place and avoid strong impact
- Please stop using this product if there is a crack or damage in it, and contact the seller or manufacturer in time

Stock condition

- Please keep it after it is deenergized
- It is forbidden to use alcohol or benzene products to clean this product after use
- Please don't keep it near the microwave oven or instrument of high voltag
- Please don't put this product near fire or in hot places
- Please don't keep it in the wet place or under the water


Machine&Accessores

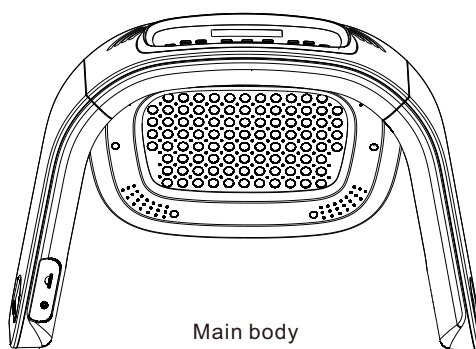Main body

Cathode bar

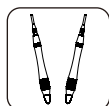
Lifting bar

Facial mask

Power adapter


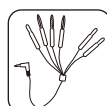
Output line


Shade goggles

Power cords are selected and matched according to standards of various countries

Parts name

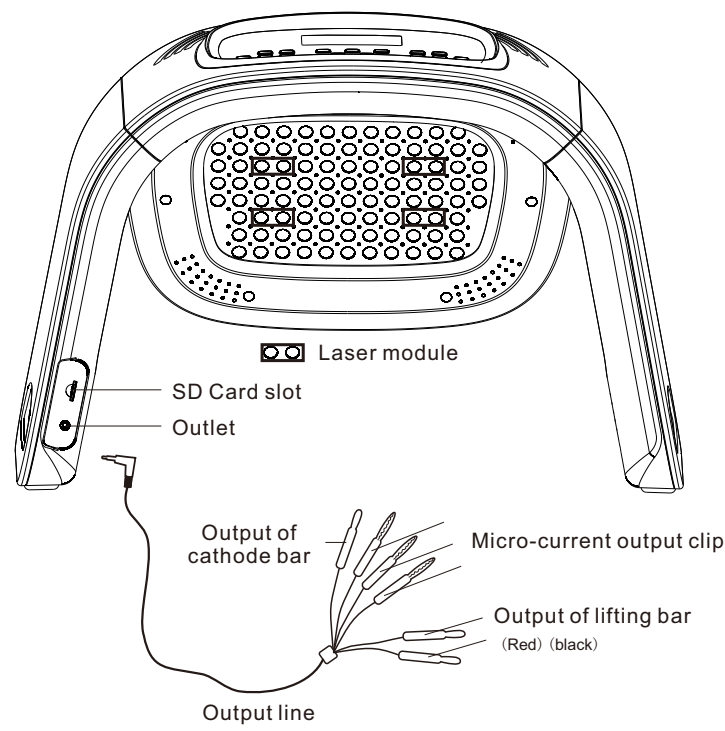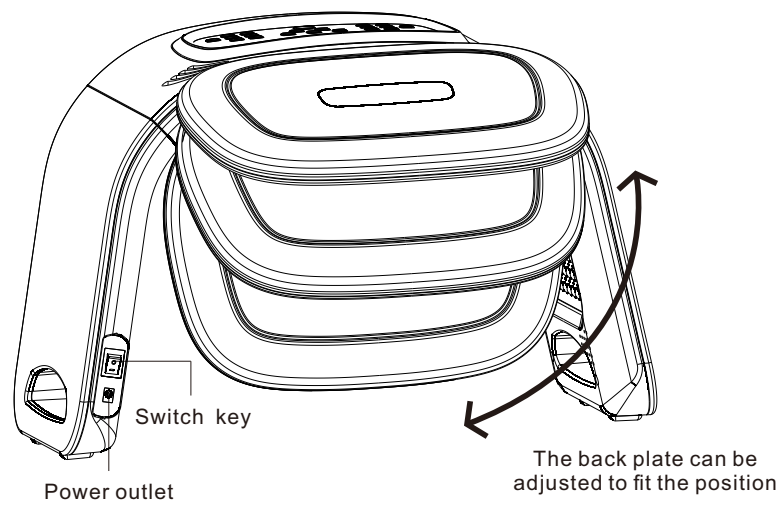Interface introduction

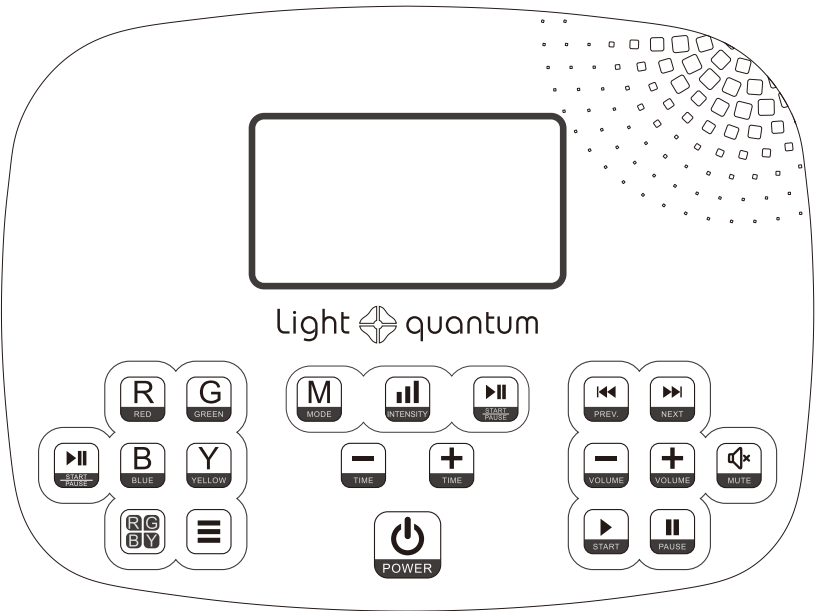
Light control	Micro-current control	Music control
Light Start button	Micro-current Mode button	Switch to the previous song
RGBY light all open	Micro-current intensity control	Switch to the next song
Light mode button	Micro-current on-off button	Mute button
Red light		Volume down
Green light		Volume up
Blue light		Start playing
Yellow light		Stop playing

Operation introduction

Neck therapy and facial therapy

Put the mask on your face and clamp the mask with the micro-current clip. Then choose the proper light color on the host interface and wear the shade patch with the cathode bar on your hand. Start the micro-current mode on the host interface and lift and pull your face with the bipolar lifting bar


(It is suggested that you should adjust from the 5th grade intensity when using micro-current mode 1-4, adjust from the 3rd grade intensity when using the micro-current mode 5, control the energy within the range of 5-20, and pull and lift your skin with the lifting bar for about 5-8 minutes)

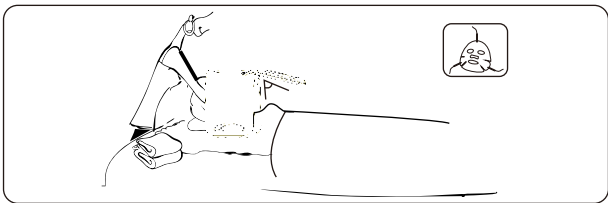It is suggested that you should lift your skin with the energy adjusted within the range of 1st -15th grade when you don't clamp the mask with the micro-current clip.

Micro-current mode: Mode 1-4 has intermittent energy output, not continuous output. It is suggested to adjust from the 3rd grade intensity.


Micro-current mode: Mode 5 has continuous output. It is suggested to adjust from the 1st grade intensity

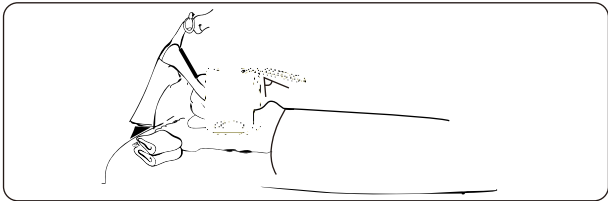Facial therapy and hair growth care

After you finish the lifting mode with gel or essence (It is suggested to use intermittent mode 1-4 and adjust from energy 2), put the mask on your face and clamp the mask with the micro-current clip. Then start the micro-current mode on the host interface with the cathode bar on your hand.

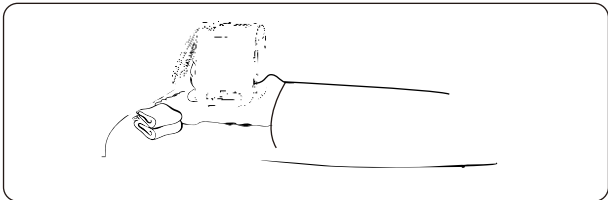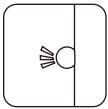
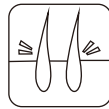
Choose the proper light color and wear the shade patch. (It is suggested that the duration of the lighting mode should be controlled in about 15-30 minutes.)

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Newface USA Inc Producer